

Certify Sulfur and Monitor Critical Elements at Sub-ppm Levels



Benchtop, Process and Portable Analysis Solutions

Sulfur Measurement Advanced

The Petra series delivers high-precision D4294 sulfur analysis across a broad measurement range. Petra MAX™ delivers D4294 sulfur analysis in addition to 12 elements from P to Zn, for rapid monitoring of critical elements like Ca, Fe, K, Ni, and V at sub-ppm levels.

Petra **MAX**
Elemental Analyzer for Petroleum

Petra **4294**
Sulfur Analyzer for Petroleum



Technical Specifications							
Petra MAX	Dynamic Range	Sulfur 5.7 ppm – 10 wt%					
	Limit of Detection (ppm @ 600 s)	Sulfur 5.7 ppm					
		P	Cl	K	Ca	V	Cr
		17	3	0.7	0.4	0.1	0.09
		Mn	Fe	Co	Ni	Cu	Zn
	0.07	0.07	0.07	0.04	0.1	0.1	
	Applications	S plus P to Zn in hydrocarbons like crude oil, diesel, gasoline, jet fuel, and lubricants.					
Petra 4294	Dynamic Range	Sulfur 2.6 ppm – 10 wt%					
	Limit of Detection (ppm @ 600 s)	Sulfur 2.6 ppm					
	Applications	S in hydrocarbons like crude oil, diesel, gasoline, jet fuel, and lubricants.					

Petra is powered by High Definition X-ray Fluorescence (HDXRF®) technology: an elemental analysis technique offering significantly enhanced detection performance over traditional XRF technology.

Advanced Workflow

Petra Series Autosampler boasts a novel design with advanced software features for a more flexible and efficient workflow. Using unique identifier (X-ID) sample cups and an open-ended sample slide, the autosampler offers sample tracking and continuous sample loading. It is an optional add-on feature for a Petra 4294 or Petra MAX analyzer. QR/barcode scanner included with purchase.



X-ID
Sample Cup

ASTM D4294
ISO 8754 | IP 336

TWO critical measurements, ONE push of a button, ZERO hassle

Sindie® +Cl is a two-in-one instrument enabling trace analysis of both sulfur and chlorine with one push of a button. It is the ideal solution to certify sulfur levels in finished products and assess chlorine for corrosion mitigation.

Features and Benefits

- Powered by MWDXRF
- LOD:
 - Sulfur: 0.4 ppm at 300 s
 - Chlorine: 0.3 ppm at 300 s
- Dynamic Range:
 - Sulfur: 0.4 ppm to 5 wt%
 - Chlorine: 0.3 ppm to 3000 ppm
- Extremely low maintenance: no conversion gasses, heating elements, columns, or quartz tubing

Options

- LIMS data output compatible software

16 32.06
Sindie + Cl
Sulfur and Chlorine Analyzer
17 35.45
Chlorine



ASTM D2622, D7039, D7536, D4929
SH / T 0842

Sulfur PLUS Lead Analysis in Petroleum Products

As aviation industry requirements become more stringent, lead detection is critical in such products as aviation gasoline and jet fuel. Sindie® +Pb delivers unprecedented precision and accuracy for sulfur and lead detection.

Features and Benefits

- Powered by MWDXRF
- LOD:
 - Sulfur: 0.7 ppm at 300 s
 - Lead: 0.002 g/L at 300 s
- Dynamic Range:
 - Sulfur: 0.7 ppm to 5000 ppm
 - Lead: 0.002 g/L to 10.000 g/L

Options

- LIMS data output compatible software

16 32.06
Sindie + Pb
Sulfur and Lead Analyzer
62 207.2
Lead



ASTM D7039 and
ISO 20884

Sulfur Analysis with Compliance Flexibility


The Sindie® 2622 series complies with ASTM D2622, D7039 and ISO 20884 methods, enabling complete flexibility in sulfur analysis. With no compromises in detection, performance and reliability, Sindie 2622 is the ideal sulfur analytical solution from ultra low sulfur diesel and gasoline to heavy fuel oil and crudes.

Features and Benefits

- Powered by MWDXRF
- Easy to use
 - Intuitive touch screen
 - Just plug-in and measure
 - Measurement time: 30-900 s
- Low and high range calibrations available:
 - Low Range: LOD - 3000 ppm
 - High Range: 0.3 wt% - 10 wt%
- Extremely low maintenance: no conversion gasses, heating elements, columns, or quartz tubing
- 75 W air-cooled excitation tube

Options

- 8-cell Autosampler
- Helium optical path available
- LIMS data output compatible software

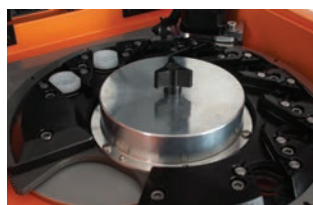

AUTOSAMPLER

- 8 sample cell capacity
- Increases productivity
- Utilizes XOS Accucell sample cups

Technical Specifications		
Sindie 2622 Gen 2	Dynamic Range	0.4 ppm to 10 wt%
	Limit of Detection (LOD)	0.4 ppm at 300 s
	Sample Cup	Traditional XRF
Sindie 2622 Gen 3	Dynamic Range	0.15 ppm to 10 wt%
	Limit of Detection (LOD)	0.15 ppm at 300 s
	Sample Cup	Accucell

ASTM D2622 & D7039
ISO 20884

Sulfur Analysis in Liquid Hydrocarbons

From ultra low sulfur diesel and gasoline, to heavy fuel oil and crudes, the Sindie® 7039 series delivers unprecedented precision and accuracy. Sindie 7039 is the ideal analytical solution for the refining industry where detection, performance and reliability are critical.

Features and Benefits

- Powered by MWDXRF
- Easy to use
 - Intuitive touch screen
 - Just plug-in and measure
 - Measurement time: 30-900 s
- Extremely low maintenance: no conversion gasses, heating elements, columns, or quartz tubing
- 75 W air-cooled excitation tube

Options

- Extended Range (XR): 0.3 wt% - 10 wt%
- LIMS data output compatible software



Technical Specifications		
Sindie 7039 Gen 2	Dynamic Range	0.4 ppm to 3000 ppm
	Limit of Detection	0.4 ppm at 300 s
	Sample Cup	Traditional XRF
Sindie 7039 Gen 3	Dynamic Range	0.15 ppm to 3000 ppm
	Limit of Detection	0.15 ppm at 300 s
	Sample Cup	Accucell
	Options	8-cell Autosampler

ASTM D7039 and
ISO 20884

Advanced Analysis with MWDXRF

Monochromatic Wavelength Dispersive X-ray Fluorescence (MWDXRF) utilizes state-of-the-art focusing and monochromating optics to increase excitation intensity and dramatically improve signal-to-background ratio compared to traditional WDXRF instruments. This enables significantly improved detection limits, precision, and a reduced sensitivity to matrix effects. A monochromatic and focused primary beam excites the sample and secondary characteristic fluorescence X-rays are emitted from the sample. A second monochromating optic selects the sulfur characteristic X-rays and directs these X-rays to the detector.

Portable Sulfur Analysis

Sindie® OTG is a portable sulfur analyzer, delivering reliable and precise results in hydrocarbons from ULSD and gasoline to marine fuels and crudes. Operate in-the-field, on-board marine environments, and in laboratories.

Features and Benefits

- Powered by MWDXRF
- LOD: 0.7 ppm at 900 s
- Dynamic Range: 0.7 ppm - 10 wt%
- Use Accucells for hassle-free sample prep
- Fits on any bench and compatible for use in mobile labs/vans
 - Dimensions: 34cm (w) x 23.5cm (d) x 30cm (h)
 - Weight: 18.5 kg
 - Utility: Standard wall power: 100-240 VAC at 2.2 A

16 32.06
Sindie OTG
Sulfur Analyzer



ASTM D7039 and
ISO 20884

Effective Online Sulfur Analysis in Petroleum Process Streams

Sindie® Online is an industrial grade process sulfur analyzer with breakthrough detection capability to monitor ultra low sulfur in petroleum or aqueous process streams. This process analyzer presents the ultimate solution for refineries and pipeline terminals where measurement speed and reliability are essential.

Features and Benefits

- Powered by MWDXRF
- Uses ASTM D7039 technology
- ATEX Zone 1 and NEC CII Div 2 Certified
- LOD: 0.5 ppmw in hydrocarbon matrices @ 300 s
- LOD: 1.5 ppmw in aqueous streams @ 300 s
- Dynamic Range: 0.5 ppmw – 3000 ppmw

Options

- Multi-stream analysis capability
- Extended Range (XR) available for measurements above 3000 ppmw up to weight percent levels

16 32.06
Sindie ONLINE
Sulfur Analyzer



ATEX and NEC Certified

XOS®
better analysis counts

15 Tech Valley Drive East Greenbush, New York 12061 USA
PH 518.880.1500 FAX 518.880.1510 info@xos.com www.xos.com

©XOS all rights reserved.
Petra MAX, Petra 4294, Sindie, HDXRF and MWDXRF are trademarks of XOS.

XOS®